

ULTIMATE



CACI NON SURGICAL FACIAL TONING & MICRODERMABRASION



ULTIMATE

Combined Non Surgical Facial Toning & Microdermabrasion System

The result of over 20 years research, the award winning CACI Facial Toning treatment has been voted* the most effective anti-ageing treatment available.

The CACI Ultimate is the only treatment system to combine CACI Non Surgical Facial Toning, Crystal Free orbital Microdermabrasion and LED Light Therapy. The system's multifunctional capabilities enable the operator to offer the most comprehensive range of non surgical treatment options available.



ORBITAL MICRODERMABRASION

This new exfoliation technology uses a hand piece with an abrasive tip that rotates at different speeds on an elliptical axis. A selection of single use, disposable abrasive treatment tips coated with a layer of silicon carbide are used to either superficially or deeply exfoliate the epidermis. This unique crystal free and vacuum free method of microdermabrasion overcomes the blockage, maintenance and downtime issues commonly associated with traditional microdermabrasion systems.

The hand piece also features the healing action of LED light therapy, which can be used simultaneously with orbital microdermabrasion or with a special massage head to stimulate the connective tissue, which improves circulation and lymph drainage.


ORBITAL MICRODERMABRASION BENEFITS

- Crystal free - no inhalation or irritation
- Vacuum free - no dragging or pulling of the skin
- Blockage free
- Disposable treatment tips for maximum hygiene
- Precision abrasion control
- Choice of 4 power levels
- A selection of 3 grades of abrasive tip



Innovative Crystal Free
Microdermabrasion
Combined With
Simultaneous LED Light
Therapy

ORBITAL MICRODERMABRASION

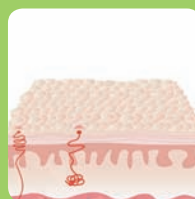

BEFORE



AFTER

The orbital rotation of the abrasive head performs a progressive and controlled skin exfoliation, removing imperfections and impurities. By removing the superficial layers of skin, the complexion is left clearer, smoother and more youthful looking. Results are visible after just one treatment.



WRINKLE COMB (LED LIGHT THERAPY)

The system's wrinkle comb combines LED photo stimulation with high frequency microcurrent. The synergy of these two technologies plumps out and softens lines, wrinkles and stretch marks, thereby providing a non invasive, needle free alternative to collagen injections and other dermal fillers.

The system has the option to select Red or Blue LED photo stimulation or an alternating pulsed LED mode that is particularly effective at treating acne and skin blemishes.




MICROCURRENT

Originally developed for treating facial palsy, microcurrent therapy is now widely used as an aesthetic anti-ageing treatment. True microcurrent uses a current with an intensity of less than one millionth of an amp (<1mA). A higher intensity current will cause contractions of the facial muscles which can worsen lines and wrinkles, whereas microcurrent stimulation does not contract muscles, but instead works by a process called muscle re-education to achieve results.

Due to its dramatic lifting and toning effect, the CACI microcurrent treatment has become known as the "Non Surgical Face Lift".

WRINKLE COMB

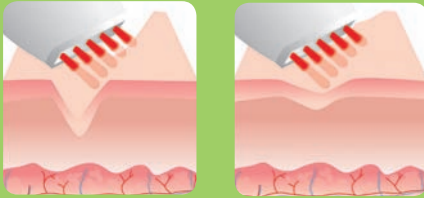

BEFORE

AFTER

Specific wavelengths of light combined with high frequency stimulation activate the body's sodium potassium pump. This encourages cell exchange, the release of sodium ions, collagen proteins and hyaluronic glycosaminoglycans which together act to plump out deep lines and wrinkles and trigger tissue repair and healing.


INCREASE YOUR REVENUE WITH CACI ULTIMATE

When you buy the CACI Ultimate, you buy more than just a treatment system, you become part of a winning team.


Share in the success of the CACI brand name and increase your revenue.

*Based on a 6 day week and 1hr treatment at £90

WEEKLY REVENUE

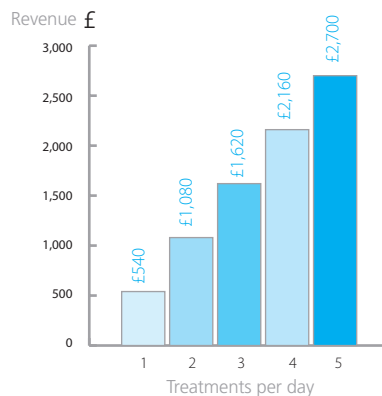

ANNUAL REVENUE

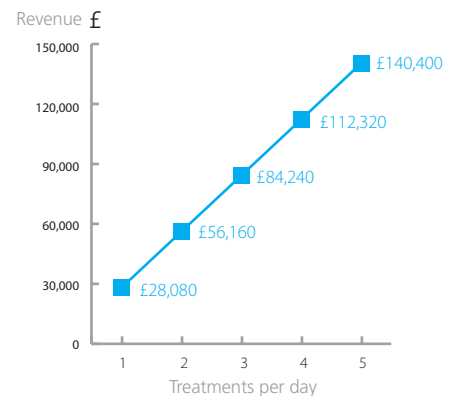



HYDRATONE

Hydratone combines active microcurrent rollers with Hydro Mask, a unique, electrically conductive silicon gel mask that has been infused with powerful hydrating properties. The rollers gently massage the face causing the mask to become electrically charged so that the whole face is bathed with rejuvenating energy. The current delivered by the rollers recharges and replenishes the bio-electrical energy within the facial muscles, restoring tone and firmness.

The energy generated also drives the hydrating actives contained within the gel mask deep into the skin tissue. The gel mask has a dramatic softening action helping to smooth lines and wrinkles. The Hydratone procedure takes just 10 minutes to perform yet provides immediate and intensive skin hydration and face firming.




86 PRESET PROGRAMS

► TREATMENTS

- Facial Toning ✓
- Lines / Wrinkles ✓
- Eye Bags ✓
- Skin Hydration ✓
- Skin Peeling ✓
- Acne / Blemishes ✓
- Neck Treatments ✓
- Cellulite (with ECM)* ✓
- Stretch Marks ✓
- Hand Treatments ✓

* Optional Attachment

► ELECTRO BUDS

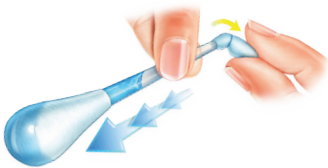

CACI systems utilise specially developed serum filled cotton buds to enhance treatment results. The serum contains Hexapeptides and powerful anti-ageing ingredients to provide an instant face firming effect. Electro buds are also perfect for retail.

ELECTRO CELLULITE MASSAGER (ECM)



The ECM is an optional attachment that can be used in conjunction with all CACI systems that incorporate microcurrent technology. The ECM attachment combines the therapeutic benefits of massage with the effectiveness of microcurrent therapy.

Conductive rollers are electrically charged with microcurrent to gently knead tired shoulders, soothe aching limbs and firm sagging muscles. The stimulatory action of the ECM tones muscles and reduces inches by improving lymph drainage and breaking down fatty deposits.

“ I am writing to express my delight and satisfaction with your staff. It is a pleasure doing business with CACI International thanks to its professional team. ”

Mrs Z Pintar
Bodyzone, Cambridge



THE CACI DIFFERENCE

All CACI systems deliver microcurrent energy using ultra low frequencies and the power of a unique digitally simulated tidal wave.

This special digital waveform optimises results and comfort by gradually increasing in intensity to overcome skin resistance (Ohm's Law). It then rises to a peak in order to deliver the full power of the microcurrent energy deep into the muscle.



Touch Screen Interface

TECHNICAL SPECIFICATION

Dimensions(cm) 45W x 42D x 15H
Weight 9Kg
Touch Screen 10.4in

Supply Voltage 61-240 volts
Consumption 110 vA
Safety Class 1 type BF in accordance with IEC 601 and IEC 601-2-10
Channels: 4 independent
Output Voltage: 30 Volts max (loaded)

MICROCURRENT
Polarity: positive, negative & alternating
Range: 20 uA to 640 uA
Frequency: 0.1Hz - 800Hz
Waveshapes: 8
Mode: Automatic Sequential

WRINKLE COMB
Output: 0.1mA - 8mA
Output Frequency: 6000Hz
LEDs 5x Blue 430nm - 40mW
LEDs 5x Red 640nm - 30mW

ORBITAL ABRASION
Rotation Speed: 0 - 4000 rpm
Motor Voltage: 6 - 15 Volts
Amperage: 0.1 - 0.5 Amps
LEDs 6x Blue 430nm - 40mW
LEDs 6x Red 640nm - 30mW

INSIDE THE BOX

ATTACHMENTS

- 1x Pair Micro Probes
- 1x Orbital Microdermabrasion Handset
- 1x Wrinkle Comb Handset
- 1x ETR Lead
- 1x Foot Pedal
- 1x Mains Lead

ACCESSORIES

- 1x T-Bar
- 1x Pair Straight Bars
- 1x Pair Hydratone Face Rollers

- 1x Pair Double Tip Probe Inserts
- 1x Abrasive Tip Remover

CONSUMABLES

- 3x Self-Adhesive Pads -Small (Set of 4)
- 1x Self-Adhesive Pads -Large (Set of 4)
- 1x Disposable Tips - Fine (Per 20)
- 1x Disposable Tips - Medium (Per 20)
- 1x Disposable Tips -Coarse (Per 20)
- 1x Dermal Massage Heads (Per 5)
- 1x Collagen Revital Serum (30g)
- 1x Essential Cleansing Milk (300ml)

- 1x Eye Wrinkle Gel (100g)
- 1x Firming Gel Concentrate (300ml)
- 1x Hydro Moisturising Cream (300ml)
- 1x Revitalising Toning Fluid (300ml)
- 1x Pack Electro-Buds (Per 100)
- 20x Hydro Masks

PROMO/MARKETING MATERIAL

- 1x A2 Poster
- 1x A5 Client Leaflets (Per 100)