

ULTRA



CACI NON SURGICAL FACIAL TONING & SKIN REJUVENATION



ULTRA

Combined Non Surgical Facial Toning & Skin Rejuvenation

The result of over 20 years research, the award winning CACI Facial Toning treatment has been voted* the most effective anti-ageing treatment available.

The CACI Ultra combines CACI Non Surgical Facial Toning with Skin Rejuvenation. The system incorporates 4 different technologies: Microcurrent, Ultrasonic Peeling, Hydratone and LED Light Therapy stimulation within one state of the art treatment device.



WRINKLE COMB (LED LIGHT THERAPY)

The Wrinkle Comb combines LED photo stimulation with high frequency electrical stimulation. The synergy of these technologies effectively plumps out and softens deep lines, wrinkles, blemishes and stretch marks. The Wrinkle Comb provides the first non invasive needle free alternative to collagen injections and other dermal fillers.

The system has the option to select Red or Blue LED photo stimulation or an alternating pulsed LED mode that is particularly effective at treating acne. The Red LED light has a healing, bio-stimulatory action on the tissue. The Blue LED light has an antibacterial and calming effect on inflamed or irritated skin.


LED LIGHT THERAPY BENEFITS

- Stimulates tissue healing
- Calming action on inflamed skin
- Effective acne treatment



Revolutionary Skin Rejuvenation Treatment

WRINKLE COMB

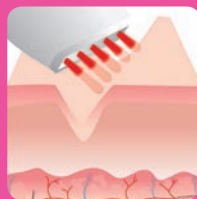

BEFORE



AFTER

Specific wavelengths of light combined with high frequency stimulation activate the body's sodium potassium pump. This encourages cell exchange, the release of sodium ions, collagen proteins and hyaluronic glycosaminoglycans which together act to plump out deep lines and wrinkles and trigger tissue repair and healing.

* Professional Beauty Magazine Survey



ULTRASONIC

Ultrasonic peeling uses a hand piece containing ceramic/quartz transducers that when excited by an electrical current emit piezoelectric vibrations. These rapid vibrations (27,000 per second) propel impurities to the skin surface and cause dead skin cells to become dislodged. The vibrating hand piece is then used to deep cleanse and gently exfoliate the epidermis. Ultrasonic skin peeling provides a far gentler and less aggressive method of skin exfoliation than traditional microdermabrasion.



HYDRATONE

Hydratone combines active microcurrent rollers with Hydro Mask, a unique, electrically conductive silicon gel mask that has been infused with powerful hydrating properties.

The rollers gently massage the face causing the mask to become electrically charged so that the whole face is bathed with rejuvenating energy. The current delivered by the rollers recharges and replenishes the bio-electrical energy within the facial muscles, restoring tone and firmness.



ULTRASONIC SKIN PEELING BENEFITS

The rapid expansions and contractions of the vibrating hand piece also destroy skin bacteria by causing the bacteria cells to implode. The Ultrasonic hand piece can also be set to a specific frequency to enhance the absorption of skin care products using sonophoresis.

- Deep cleanse & gently exfoliate the epidermis
- Destroy skin bacteria
- Enhance absorption of skin care formulations
- Stimulate and improve circulation


INCREASE YOUR REVENUE WITH CACI ULTRA

When you buy the CACI Ultra, you buy more than just a treatment system, you become part of a winning team.


Share in the success of the CACI brand name and increase your revenue.

*Based on a 6 day week and 1hr treatment at £70

WEEKLY REVENUE

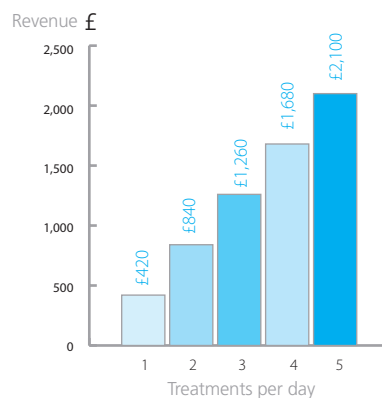

ANNUAL REVENUE

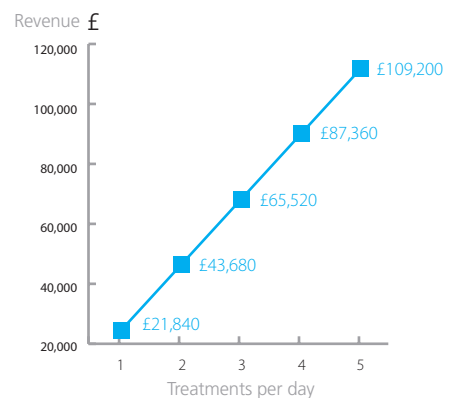



MICROCURRENT

Originally developed for treating facial palsy, microcurrent therapy is now widely used as an aesthetic anti-ageing treatment.

True microcurrent uses a current with an intensity of less than one millionth of an amp (<1mA) that works in harmony with the body's own bio-electrical field. A higher intensity current will cause contractions of the facial muscles which can actually worsen lines and wrinkles, whereas microcurrent stimulation does not contract muscles, but instead works by a process called muscle re-education to achieve results.

Due to its dramatic lifting and toning effect, the CACI microcurrent treatment has become known as the "Non Surgical Face Lift".




82 PRESET PROGRAMS

▶ TREATMENTS

Facial Toning	✓
Lines / Wrinkles	✓
Eye Bags	✓
Skin Hydration	✓
Skin Peeling	✓
Acne / Blemishes	✓
Neck Treatments	✓
Cellulite (with ECM)*	✓
Stretch Marks	✓
Hand Treatments	✓

* Optional Attachment

▶ ELECTRO BUDS

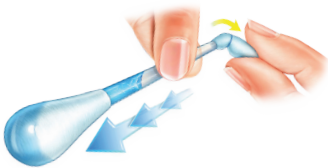

CACI systems utilise specially developed serum filled cotton buds to enhance treatment results. The serum contains Hexapeptides and powerful anti-ageing ingredients to provide an instant face firming effect. Electro buds are also perfect for retail.

ELECTRO CELLULITE MASSAGER (ECM)




The ECM is an optional attachment that can be used in conjunction with all CACI systems that incorporate microcurrent technology. The ECM attachment combines the therapeutic benefits of massage with the effectiveness of microcurrent therapy.

Conductive rollers are electrically charged with microcurrent to gently knead tired shoulders, soothe aching limbs and firm sagging muscles. The stimulatory action of the ECM tones muscles and reduces inches by improving lymph drainage and breaking down fatty deposits.

“ I have been thrilled with my Ultra machine as have my clients, we had a launch evening and I was able to pay off 50% of the outlay of the machine in this first evening. Since then we have been inundated with clients wanting treatment ”

Candice White Beauty Salon, Bristol

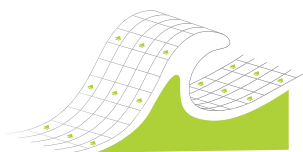

THE CACI DIFFERENCE

All CACI systems deliver microcurrent energy using ultra low frequencies and the power of a unique digitally simulated tidal wave.

This special digital waveform optimises results and comfort by gradually increasing in intensity to overcome skin resistance (Ohm's Law). It then rises to a peak in order to deliver the full power of the microcurrent energy deep into the muscle.



Touch Screen Interface

TECHNICAL SPECIFICATION

Dimensions(cm) 45W x 42D x 15H
Weight 9Kg
Touch Screen 10.4in
Supply Voltage 61-240 volts
Consumption 110 vA
Safety Class 1 type BF in accordance with IEC 601 and IEC 601-2-10

MICROCURRENT
Polarity: positive, negative & alternating
Range: 20 uA to 640 uA
Frequency: 0.1Hz - 800Hz
Waveshapes: 8
Mode: Automatic Sequential

ULTRASONIC HANDSET
Mechanic Oscillation: 30micro +/- 20%
Frequency: 27,000Hz +/- 10%

WRINKLE COMB
Output: 0.1mA - 8mA
Output Frequency: 6000Hz
LEDs 5x Blue 430nm - 40mW
LEDs 5x Red 640nm - 30mW

INSIDE THE BOX

ATTACHMENTS
 1x Pair Micro Probes
 1x Ultrasonic Handset
 1x Wrinkle Comb Handset
 1x ETR Lead
 1x Foot Pedal
 1x Mains Lead

ACCESSORIES
 1x T-Bar
 1x Pair Straight Bars
 1x Pair Hydratone Face Rollers
 1x Pair Double Tip Probe Inserts

CONSUMABLES
 3x Self-Adhesive Pads -Small (Set of 4)
 1x Self-Adhesive Pads -Large (Set of 4)
 1x Ultra Pre-Peeling Solution (250ml)
 1x Ultra Moisturising Gel (250ml)
 1x Ultra Healing Gel (250ml)
 1x Ultra Sterilising Solution (250ml)
 1x Collagen Revital Serum (30g)
 1x Essential Cleansing Milk (300ml)
 1x Eye Wrinkle Gel (100g)
 1x Firming Gel Concentrate (300ml)
 1x Hydro Moisturising Cream (300ml)
 1x Revitalising Toning Fluid (300ml)

1x Pack Electro-Buds (Per 100)
 20x Hydro Masks

PROMO/MARKETING MATERIAL
 1x A2 Poster
 1x A5 Client Leaflets (Per 100)